

# HIGHLY ACCURATE HIFI READS FOR AAV RESEARCH

PacBio® long-read HiFi sequencing allows you to sequence the entire AAV genome with extraordinary accuracy and completeness so that you can quickly discover and optimize gene therapy vectors.

**Your moment of discovery in AAV research is waiting.**

## What can you do with highly accurate AAV sequencing?

- Rapidly accelerate the engineering and discovery of novel capsids with desired characteristics such as tissue specificity through full-length AAV genome or cap gene sequencing.
- Monitor and improve AAV manufacturing by characterizing AAV vectors and their impurities with high accuracy.
  - Assess packaged impurity sequence profiles like identity, size distribution, and relative abundance.
  - Profile payload sequence length and integrity.
  - Characterize ITR rearrangements and integrity.
- Confidently evaluate and optimize your AAV design.
  - Identify payload sequence truncation hotspots.
  - Verify whether a target construct or gene of interest is correctly expressed and spliced.
- Explore potential risks by assessing host genome integration events.

PacBio highly accurate long-read sequencing offers advantages over short reads with these features and is recommended for rAAV characterization.<sup>1</sup>

- Long reads up to 25 kb span the entire genes, regions of interest, or AAV genomes, including inverted terminal repeats (ITRs).
- Q30 accuracy (99.9%) provides Sanger-quality base-level resolution.
- Detect variants beyond small insertions and deletions to ensure both base-level modifications and structural rearrangements are captured.
- Uniform coverage to detect variants in repetitive and extreme GC-rich regions in an unbiased manner.
- Sequence transcripts and evaluate splicing events with full-length cDNA sequencing.
- Easy-to-implement protocol for AAV sequencing with end-to-end workflow suitable for both scAAV and ssAAV.

## Workflow



### Sample prep

Use a range of starting materials that fit your project:

- Single-stranded ssAAV
- Self-complementary scAAV
- Target capture-enriched DNA
- DNA amplicons
- RNA



### Library prep

Use standardized and automatable workflows with PacBio Compatible partners.

Supports a range of target types and multiplexing options for up to 96 targets



### HiFi sequencing

Generate 60 Gb on Vega or up to 480 Gb per run on Revio with SPRQ™ chemistry.



### Data analysis

Use SMRT® Link, a web-based end-to-end workflow to demultiplex, analyze, and visualize your sequencing data.

Further analyze your data using publicly available bioinformatic tools,<sup>2</sup> or the Form Bio<sup>3</sup> platform.

## How can HiFi sequencing empower your AAV research?

### Help identify novel AAV capsids

Identifying novel capsids with tissue-specific tropism is critical for targeted delivery of AAV-based gene therapies.<sup>4,5</sup>

In one study, HiFi sequencing of intact full-length AAV genomes led to the discovery of novel capsid variants.<sup>6</sup>

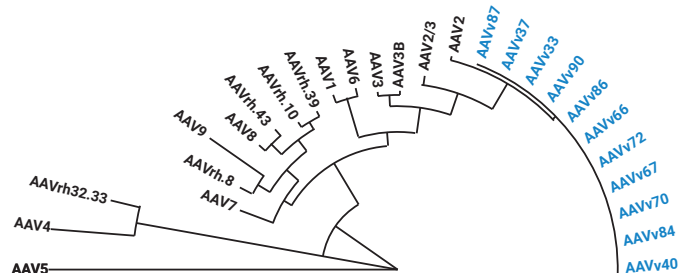


Figure 1. (from Hsu et al. 2020) Phylogenetic tree of AAV2 variants reported in this study (blue) and contemporary serotypes.

### Assess gene transfer vectors

Performing highly accurate long-read sequencing enables a more complete characterization of the viral genome to ensure correct packaging.<sup>7,8</sup>

HiFi sequencing provides more complete viral sequences for heterogeneous vector populations. This allows for unambiguous detection and classification of plasmid and host genomic encapsidated DNA.

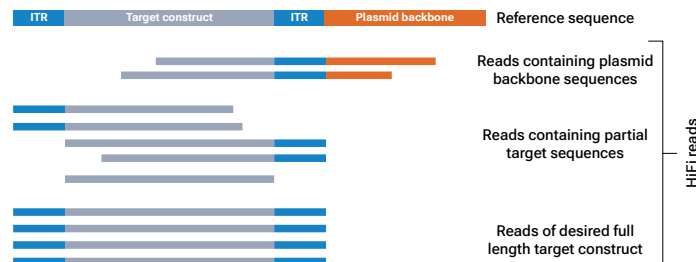


Figure 2. HiFi long reads span the full length of a viral genome including ITRs, target sequences, and plasmid backbone sequences.

### Characterize ITR integrity and configuration

ITR integrity is correlated to vector genome heterogeneity. HiFi sequencing allows for full ITR sequencing and resolution of configurations and mutations that can inform understanding of AAV genome heterogeneity.<sup>7</sup>

### Identify payload truncation hotspots

Secondary structure elements in AAV payloads such as in shRNA sequences (small hairpin RNAs) can lead to rAAV truncation hotspots.<sup>9</sup>

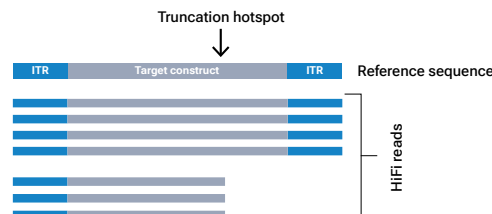


Figure 3. HiFi long reads reveal truncation hotspots within AAV payloads.

### Ensure a transcript is correctly expressed and spliced using the Iso-Seq® method

Accurate splicing can be a critical factor to consider in gene therapy research. The Iso-Seq method can capture full-length transcripts and reveal isoforms not distinguishable by short-read sequencing.<sup>10</sup>

### Help enable the assessment of AAV integration into the host cell genome, which is critical for evaluating the safety profile of the AAV vector.

HiFi sequencing allows unambiguous identification of integration sites and resolves integrated concatemers and rearrangements that cannot be detected with short-read methods.<sup>11</sup>

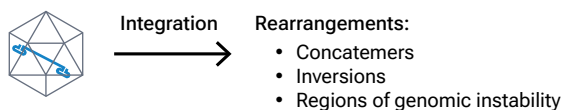


Figure 4. HiFi sequencing allows for the assessment of potential genomic rearrangements due to AAV integration.



Learn about highly accurate HiFi reads for AAV gene therapy research: [pacb.com/gene-therapy](https://pacb.com/gene-therapy)

1. Fuentes, Staudhammer J, Wright JF, et al. (2023) DHC white paper: Beyond empty and full: Understanding heterogeneity in rAAV products and impurities. <https://darkhorseconsultinggroup.com/heterogeneity-in-raav-vector-genomes-etc>
2. <https://github.com/Magdoll/AAV>
3. <https://info.formbio.com/pacbio>
4. Paulk NK, Pekrun K, Zhu E, et al (2018) Bioengineered AAV capsids with combined high human liver transduction in vivo and unique humoral seroreactivity. *Mol Ther* 26:289–303. <https://doi.org/10.1016/j.ymthe.2017.09.021>
5. Casy W, Garza IT, Chen X, et al (2023) SMRT sequencing enables high-throughput identification of novel AAVs from capsid shuffling and directed evolution. *Genes* 14:1660. <https://doi.org/10.3390/genes14081660>
6. Hsu H-L, Brown A, Loveland AB, et al (2020) Structural characterization of a novel human adeno-associated virus capsid with neurotropic properties. *Nat Commun* 11:3279. <https://doi.org/10.1038/s41467-020-17047-1>
7. Tai PWL, Xie J, Fong K, et al (2018) Adeno-associated virus genome population sequencing achieves full vector genome resolution and reveals human-vector chimeras. *Mol Ther Methods Clin Dev* 9:130–141. <https://doi.org/10.1016/j.omtm.2018.02.002>
8. Tran NT, Lecomte E, Saleun S, et al (2022) Human and insect cell-produced recombinant adeno-associated viruses show differences in genome heterogeneity. *Hum Gene Ther* 33:371–388. <https://doi.org/10.1089/hum.2022.050>
9. Xie J, Mao Q, Tai PWL, et al (2017) Short DNA hairpins compromise recombinant adeno-associated virus genome homogeneity. *Mol Ther* 25:1363–1374. <https://doi.org/10.1016/j.ymthe.2017.03.028>
10. Cripe TP, Hutzen B, Currier MA, et al (2022) Leveraging gene therapy to achieve long-term continuous or controllable expression of biotherapeutics. *Sci Adv* 8:eabm1890. <https://doi.org/10.1126/sciadv.abm1890>
11. Dalwadi DA, Calabria A, Tiyaboonchai A, et al (2021) AAV integration in human hepatocytes. *Mol Ther* 29:2898–2909. <https://doi.org/10.1016/j.ymthe.2021.08.031>

Research use only. Not for use in diagnostic procedures. © 2025 Pacific Biosciences of California, Inc. ("PacBio"). All rights reserved. Information in this document is subject to change without notice. PacBio assumes no responsibility for any errors or omissions in this document. Certain notices, terms, conditions and/or use restrictions may pertain to your use of PacBio products and/or third-party products. Refer to the applicable PacBio terms and conditions of sale and to the applicable license terms at [pacb.com/license](https://pacb.com/license). Pacific Biosciences, the PacBio logo, Circulomics, Omniome, SMRT, SMRTbell, Iso-Seq, Sequel, Nanobind, SBB, Revio, Onso, Apton, Kinnex, PureTarget, SPRQ, and Vega are trademarks of PacBio.