

# Procedure & Checklist – cDNA Capture Using IDT xGen® Lockdown® Probes

## Before You Begin


This document describes the process for capturing cDNA prepared with the SMARTer PCR cDNA Synthesis Kit (Clontech) and pulled-down using xGen Lockdown Probes/Panels from IDT.

To perform this procedure, you must have reviewed the [Procedure & Checklist – Iso-Seq Template Preparation for Sequel Systems](#).

## Workflow

The workflow includes the following:

1. Preparing the cDNA library using the SMARTer PCR cDNA Synthesis Kit.
2. Capturing cDNA with the IDT xGen Lockdown Probes (biotinylated probes).
3. Constructing SMRTbell® libraries
4. Sequencing using the PacBio System.

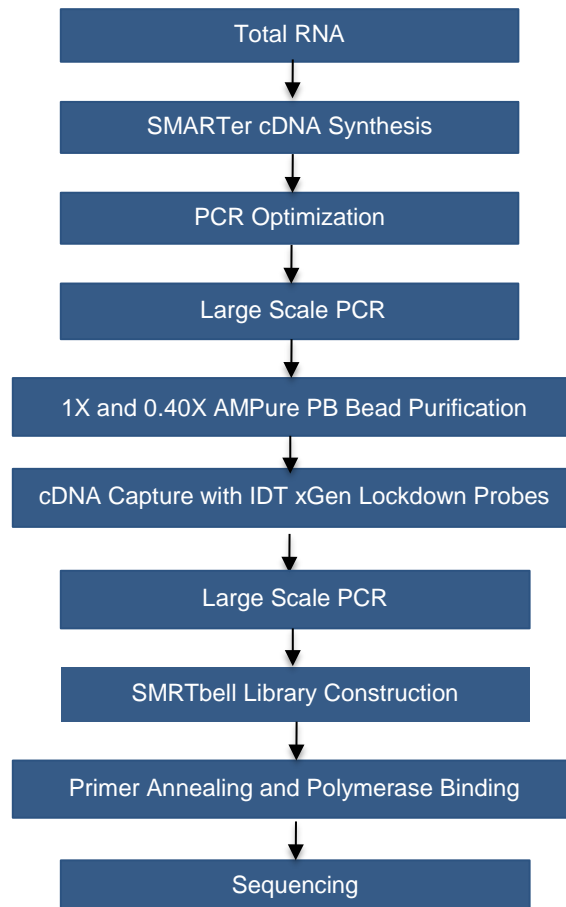

## Materials Needed

Item	Vendor	Part Number
<b>cDNA Library</b> SMARTer PCR cDNA Synthesis Kit Lonza flash gel	KAPA Biosystems Lonza	KK8503
<b>PCR/Target Capture</b> PolyT blocker Oligo (5' TTT TTT TTT TTT TTT TTT TTT TTT TTT / 3InvdT/3') xGen Hybridization and Wash Kit (16 or 96 reaction) xGen Lockdown Panels/Probes (target probes)* Takara LA Taq DNA Polymerase Hot-Start version SMARTer PCR Oligo (5' AAG CAG TGG TAT CAA CGC AGA GTA C 3') Dynabeads® M-270 Streptavidin AMPure® PB beads	IDT IDT IDT Clontech IDT Life Technologies PacBio	N/A N/A RR042A N/A 65305
<b>SMRTbell Library Construction and Sequencing</b> Template Prep Kit DNA/Polymerase Binding Kit DNA Sequencing Kit AMPure® PB beads	PacBio PacBio PacBio	

\*If the xGen Lockdown Probes were received dry, resuspend them in IDTE pH 8.0 to a final concentration of 0.75pmol/μL. For additional support regarding resuspension of Lockdown Probes pools, visit [www.idtdna.com/xgen](http://www.idtdna.com/xgen) ► xGen Lockdown Probes ► Support tab ► expand Number of Reactions and Resuspension Volume.

## Prepare cDNA Library

To prepare a cDNA library, refer to pages 1 – 11 ([Procedure & Checklist – Iso-Seq Template Preparation for Sequel Systems](#)).

1. Prepare 1st-strand synthesis using the SMARTer PCR cDNA Synthesis Kit.
2. Enrich by:
  - a. Optimizing PCR cycles.
  - b. Performing large-scale PCR.

STEP	✓	Prepare the Hybridization Sample	Notes								
		<p>In this section, you will need the following:</p> <ul style="list-style-type: none"> <li>• SMARTer PCR oligo (IDT)</li> <li>• PolyT blocker (IDT)</li> <li>• 2X Hybridization Buffer contained in xGen Lockdown Hybridization and Wash Kit</li> <li>• Hybridization Buffer Enhancer in xGen Lockdown Hybridization and Wash Kit</li> <li>• xGen Lockdown Panels/Probes (target probes)</li> </ul>									
<b>1</b>		Add 1 – 1.5 µg cDNA to a new 0.5 mL LoBind tube.									
<b>2</b>		Add 1 µL of SMARTer PCR oligo and 1 µL PolyT blocker (both at 1000 µM) to the tube containing the cDNA.									
<b>3</b>		Close the tube's lid and puncture a hole in the cap with an 18 – 20 gauge or smaller needle.									
<b>4</b>		Dry the cDNA Sample Library/SMARTer PCR oligo/PolyT blocker completely in a LoBind tube using a DNA vacuum concentrator (speed vac). Do not leave tubes in the speed vac once they have dried. This will result in over drying the tube contents.									
<b>5</b>		<p>To the dried-down sample add:</p> <table border="1" data-bbox="448 932 1179 1131"> <thead> <tr> <th>Component</th> <th>Volume</th> </tr> </thead> <tbody> <tr> <td>2X Hybridization Buffer</td> <td>8.5 µL</td> </tr> <tr> <td>Hybridization Buffer Enhancer</td> <td>2.7 µL</td> </tr> <tr> <td>Nuclease Free Water</td> <td>1.8 µL</td> </tr> </tbody> </table>	Component	Volume	2X Hybridization Buffer	8.5 µL	Hybridization Buffer Enhancer	2.7 µL	Nuclease Free Water	1.8 µL	
Component	Volume										
2X Hybridization Buffer	8.5 µL										
Hybridization Buffer Enhancer	2.7 µL										
Nuclease Free Water	1.8 µL										
<b>6</b>		Seal the hole in the tube's cap with a laboratory tape.									
<b>7</b>		Mix the reaction by tapping the tube, followed by a quick spin.									
<b>8</b>		Place the tube in a +95°C heat block for 10 minutes to denature the cDNA.									
<b>9</b>		Quick spin at maximum speed, allowing the mix to cool to room temperature before addition of probes. Probes should never be added while at 95°C.									
<b>10</b>		Leave the PCR tube on the bench for approximately 30 seconds, then add 4 µL of xGen Lockdown Panel/Probe for a total volume of 17 µL.									
<b>11</b>		Mix and quick spin.									
<b>12</b>		Incubate in a thermocycler at +65°C for 4 hours. The thermocycler's heated lid should be turned on and set to maintain +75°C (10°C above the hybridization temperature).									

STEP	✓	Preparing Beads for Capture	Notes																																				
		<p>In this section, you will need the following:</p> <ul style="list-style-type: none"> <li>Wash buffers contained in the xGen Lockdown Hybridization and Wash Kit</li> <li>Dynabeads M-270 Streptavidin</li> </ul>																																					
1		<p>Prepare Wash Buffers:</p> <p>a. Prepare 1X working solutions of the buffers listed in the below table. The total volume of 1X buffer in the table is for a single experiment. Scale up accordingly when multiple samples are required.</p> <p>Ensure that 10x Wash Buffer I is in solution before use. Any solids can be dissolved by warming with your hands or incubate at 65°C heat block and vortexing.</p> <table border="1" data-bbox="423 537 1341 848"> <thead> <tr> <th>Buffer Stock</th> <th>Stock Conc.</th> <th>Vol. Buffer</th> <th>Vol. Water</th> <th>Total Volume*</th> <th>Final Conc.</th> </tr> </thead> <tbody> <tr> <td>Wash Buffer I (tube 1)</td> <td>10X</td> <td>30 µL</td> <td>270 µL</td> <td>300 µL</td> <td>1X</td> </tr> <tr> <td>Wash Buffer II (tube 2)</td> <td>10X</td> <td>20 µL</td> <td>180 µL</td> <td>200 µL</td> <td>1X</td> </tr> <tr> <td>Wash Buffer III (tube 3)</td> <td>10X</td> <td>20 µL</td> <td>180 µL</td> <td>200 µL</td> <td>1X</td> </tr> <tr> <td>Stringent Wash Buffer (tube S)</td> <td>10X</td> <td>40 µL</td> <td>360 µL</td> <td>400 µL</td> <td>1X</td> </tr> <tr> <td>Bead Wash Buffer</td> <td>2X</td> <td>250 µL</td> <td>250 µL</td> <td>500 µL</td> <td>1X</td> </tr> </tbody> </table> <p>*Store working solutions at room temperature (+15 to +25°C) for up to 4 weeks. The volumes in this table are calculated for a single capture; scale up accordingly if multiple hybridization reactions will be processed.</p> <p>b. Preheat the following wash buffers to +65°C in a heat block or water bath:</p> <ul style="list-style-type: none"> <li>100 µL of 1X Wash Buffer I (Tube 1)</li> <li>400 µL of 1X Stringent Wash Buffer (Tube S)</li> </ul>	Buffer Stock	Stock Conc.	Vol. Buffer	Vol. Water	Total Volume*	Final Conc.	Wash Buffer I (tube 1)	10X	30 µL	270 µL	300 µL	1X	Wash Buffer II (tube 2)	10X	20 µL	180 µL	200 µL	1X	Wash Buffer III (tube 3)	10X	20 µL	180 µL	200 µL	1X	Stringent Wash Buffer (tube S)	10X	40 µL	360 µL	400 µL	1X	Bead Wash Buffer	2X	250 µL	250 µL	500 µL	1X	
Buffer Stock	Stock Conc.	Vol. Buffer	Vol. Water	Total Volume*	Final Conc.																																		
Wash Buffer I (tube 1)	10X	30 µL	270 µL	300 µL	1X																																		
Wash Buffer II (tube 2)	10X	20 µL	180 µL	200 µL	1X																																		
Wash Buffer III (tube 3)	10X	20 µL	180 µL	200 µL	1X																																		
Stringent Wash Buffer (tube S)	10X	40 µL	360 µL	400 µL	1X																																		
Bead Wash Buffer	2X	250 µL	250 µL	500 µL	1X																																		
2		<p>Prepare the capture beads:</p> <ol style="list-style-type: none"> <li>Allow the Dynabeads M-270 Streptavidin to warm to room temperature for 30 minutes prior to use.</li> <li>Mix the beads thoroughly by vortexing for 15 seconds.</li> <li>For a single sample, aliquot 100 µL beads into a 1.5 mL LoBind tube. Scale up volume for multiple samples.</li> <li>Place the LoBind tube in a magnetic rack. When the supernatant is clear, remove and discard the supernatant being careful not to disturb the beads. Any remaining traces of liquid will be removed with subsequent wash steps. Note: Allow the Dynabeads to settle for at least 1-2 minutes before removing the supernatant. The Dynabeads are “filmy” and slow to collect to the side of the tube.</li> <li>While the LoBind tube is in the magnetic rack, add 200 µL of 1X Bead Wash Buffer. For multiple samples, wash with 200 µL x X samples.</li> <li>Remove the tube from the magnetic rack and vortex until the beads are in solution.</li> <li>Quickly spin and place the LoBind tube back in the magnetic rack to collect the beads to the side of the tube. Once clear, remove and discard the liquid.</li> <li>Repeat steps e - g for a total of two washes.</li> <li>Resuspend by vortexing the beads in 100 µL of 1X Bead Wash Buffer. For multiple samples, scale up accordingly.</li> <li>Place the tube in the magnetic rack to collect beads to the side of the tube. For multiple samples, transfer 100 µL aliquot into new LoBind tube. Once clear, remove and discard the supernatant.</li> <li>The washed beads are now ready to bind the captured DNA. Proceed immediately to the next step. Do not allow the capture beads to dry. Small amounts of residual Bead Wash Buffer will not interfere with binding of DNA to the capture beads.</li> </ol>																																					

STEP	Binding cDNA to Beads and Wash	Notes
1	Bind cDNA to the capture beads: <ol style="list-style-type: none"> <li>a. To the washed capture beads, transfer the 17 <math>\mu</math>L hybridized probe/sample mixture prepared in the “Preparing hybridization section”.</li> <li>b. Mix by tapping the tube until the sample is homogeneous.</li> <li>c. Incubate in a thermomixer set to +65°C for 45 minutes or transfer the mix to a PCR tube and incubate in a thermocycler (heated lid set to +75°C). Hand mix periodically by gently tapping the tube to keep the beads in suspension.</li> </ol>	
2	Wash the captured cDNA: <ol style="list-style-type: none"> <li>a. Pre-heat 1X Wash Buffer (tube 1) and 1X Stringent Wash Buffer (Tube S) to 65°C,</li> <li>b. After 45 minutes of incubation, remove the tube from the 65°C thermomixer and add 100 <math>\mu</math>L pre-heated 1X Wash Buffer I (Tube 1).</li> <li>c. Mix thoroughly by tapping the tube until the sample is homogeneous.</li> <li>d. If using a PCR tube, transfer the sample to a 1.5 mL LoBind tube.</li> <li>e. Place the tube in the magnetic rack to collect the beads to the side of the tube. Remove and discard the liquid once clear.</li> <li>f. Remove the tube from the magnetic rack and add 200 <math>\mu</math>L of 1X Stringent Wash Buffer (TubeS) heated to +65°C. Mix by tapping the tube until the sample is homogeneous. Work quickly so that the temperature does not drop below +65°C.</li> <li>g. Incubate at +65°C for 5 minutes.</li> <li>h. Repeat steps d - f for a total of two washes using 1X Stringent Wash Buffer (TubeS) heated to +65°C.</li> <li>i. Place the tubes in the magnetic rack to collect the beads to the side of the tube. Remove and discard the liquid once clear.</li> <li>j. Add 200 <math>\mu</math>L of room temperature 1X Wash Buffer I (Tube1). Hand mix by gently tapping the tube. Quick spin.</li> <li>k. Place the tube in the magnetic rack to collect the beads to the side of the tube. Remove and discard the liquid once clear.</li> <li>l. Add 200 <math>\mu</math>L of room temperature 1X Wash Buffer II (Tube2) and mix thoroughly by tapping the tube until Sample is homogeneous. Quick Spin.</li> <li>m. Place the tube in the magnetic rack to collect the beads to the side of the tube. Remove and discard the liquid once clear.</li> <li>n. Add 200 <math>\mu</math>L of room temperature 1X Wash Buffer III (Tube 3) and mix by tapping the tube until sample is homogeneous. Quick Spin.</li> <li>o. Place the tubes in the magnetic rack to collect the beads to the side of the tube. Remove and discard the liquid once clear.</li> <li>p. Remove the tubes from the magnetic rack and add 50 <math>\mu</math>L of EB to each tube of bead-bound captured sample. This is enough for two PCR reactions required in the next section.</li> <li>q. Store the beads plus captured samples at -15 to -25°C or proceed to the next step. It is not necessary to separate the beads from the eluted DNA. The bead/sample mix can be added to the PCR reaction directly.</li> </ol>	






STEP	✓	Amplification of Captured DNA Sample	Notes																																					
		In this section, you will need the following: <ul style="list-style-type: none"> <li>• Takara LA Taq DNA Polymerase Hot-Start Version from Clontech</li> <li>• SMARTer PCR Oligos (from the Clontech Kit)</li> </ul>																																						
1		PacBio recommends using Takara LA Taq DNA Polymerase Hot-Start version. <p>a. Assemble the following PCR reaction:</p> <table border="1" data-bbox="461 499 1136 884"> <thead> <tr> <th>Component</th> <th>Volume</th> </tr> </thead> <tbody> <tr> <td>Water</td> <td>104.5 µL</td> </tr> <tr> <td>10x LA PCR Buffer</td> <td>20 µL</td> </tr> <tr> <td>2.5 mM each dNTPs</td> <td>16 µL</td> </tr> <tr> <td>SMARTer PCR Oligos (12 µM each)</td> <td>8.3 µL</td> </tr> <tr> <td>Takara LA Taq DNA polymerase</td> <td>1.2 µL</td> </tr> <tr> <td>Captured Library</td> <td>50 µL</td> </tr> <tr> <td>Total</td> <td>200 µL</td> </tr> </tbody> </table> <p>b. Split the PCR mix into two tubes, 100 µL each. It is best to perform the PCR reaction in 100 µL volumes.</p> <p>c. Amplify using the following PCR protocol:</p> <table border="1" data-bbox="461 999 1143 1362"> <thead> <tr> <th>Step</th> <th>Temp</th> <th>Time</th> </tr> </thead> <tbody> <tr> <td>1</td> <td>95°C</td> <td>2 minutes</td> </tr> <tr> <td>2</td> <td>95°C</td> <td>20 seconds</td> </tr> <tr> <td>3</td> <td>68°C</td> <td>10 minutes</td> </tr> <tr> <td>4</td> <td colspan="2">Repeat steps 2-3, 7 to 10 times for a total of 8 to 11 cycles</td> </tr> <tr> <td>5</td> <td>72°C</td> <td>10 minutes</td> </tr> <tr> <td>6</td> <td>4°C</td> <td>Hold</td> </tr> </tbody> </table>	Component	Volume	Water	104.5 µL	10x LA PCR Buffer	20 µL	2.5 mM each dNTPs	16 µL	SMARTer PCR Oligos (12 µM each)	8.3 µL	Takara LA Taq DNA polymerase	1.2 µL	Captured Library	50 µL	Total	200 µL	Step	Temp	Time	1	95°C	2 minutes	2	95°C	20 seconds	3	68°C	10 minutes	4	Repeat steps 2-3, 7 to 10 times for a total of 8 to 11 cycles		5	72°C	10 minutes	6	4°C	Hold	
Component	Volume																																							
Water	104.5 µL																																							
10x LA PCR Buffer	20 µL																																							
2.5 mM each dNTPs	16 µL																																							
SMARTer PCR Oligos (12 µM each)	8.3 µL																																							
Takara LA Taq DNA polymerase	1.2 µL																																							
Captured Library	50 µL																																							
Total	200 µL																																							
Step	Temp	Time																																						
1	95°C	2 minutes																																						
2	95°C	20 seconds																																						
3	68°C	10 minutes																																						
4	Repeat steps 2-3, 7 to 10 times for a total of 8 to 11 cycles																																							
5	72°C	10 minutes																																						
6	4°C	Hold																																						
2		After amplification, pool the 100 µL reactions and proceed to the next step to purify products using AMPure PB beads.																																						

STEP	✓	Post Amplification Clean Up	Notes
1		Add <b>1X</b> AMPure PB beads to the pooled PCR product.	
2		Mix by tapping the LoBind tube until the sample is homogeneous.	
3		Incubate at room temperature for 10 minutes.	
4		Place on magnetic rack until solution clears. Remove and discard supernatant.	
6		With the tube still on magnet, add 200 $\mu$ L freshly prepared 70% ethanol to the tube containing beads plus DNA.	
7		Remove and discard 70% ethanol.	
8		Repeat steps 5 to 6 for total of two washes with 70% ethanol.	
9		Let beads air dry for 1 minute. (Note - over drying the beads will result in reduced DNA yield.)	
10		Add 27 $\mu$ L EB and remove the tube from the magnet. Mix by tapping the tube until the sample is homogeneous. Then incubate at room temperature for 2 minutes.	
11		Place back on magnet. When the solution clears, remove 25 $\mu$ L supernatant into new 1.5 mL LoBind tube.	
12		Determine concentration using Qubit device or similar quantification assay.	
13		Run 1 $\mu$ L of sample on Agilent DNA 12000 chip according to manufacturer's instructions.	
14		The captured cDNA is now ready for SMRTbell library construction.	

## Repair DNA Damage

Use the following table to repair any DNA damage. If preparing larger amounts of DNA, scale the reaction volumes accordingly (i.e., for 10 µg of DNA scale the total volume to 100 µL). Do not exceed 100 ng/µL of DNA in the final reaction.

1. In a LoBind microcentrifuge tube, add the following reagents:


Reagent	Cap Color	Stock Conc.	Volume	Final Conc.	✓	Notes
cDNA			__ µL for 5.0 µg			
DNA Damage Repair Buffer		10 X	5.0 µL	1 X		
NAD <sup>+</sup>		100 X	0.5 µL	1 X		
ATP high		10 mM	5.0 µL	1 mM		
dNTP		10 mM	0.5 µL	0.1 mM		
DNA Damage Repair Mix			2.0 µL			
H <sub>2</sub> O			__ µL to adjust to 50.0* µL			
<b>Total Volume</b>			50.0 µL			

\*To determine the correct amount of H<sub>2</sub>O to add, use your actual DNA amount noted in the Notes column.

2. Mix the reaction well by gentle mixing.
3. Spin down contents of LoBind tube with a quick spin in a microfuge.
4. Incubate at 37°C for 20 minutes, then return the reaction to 4°C for 1 minute.

## Repair Ends

Use the following table to prepare your reaction then purify the DNA.

Reagent	Tube Cap Color	Stock Conc.	Volume	Final Conc.	✓	Notes
cDNA (Damage Repaired)			50 µL			
End Repair Mix		20 X	2.5 µL	1X		
<b>Total Volume</b>			52.5 µL			

1. Mix the reaction well by gentle mixing.
2. Spin down contents of LoBind tube with a quick spin in a microfuge.
3. Incubate at 25°C for 5 minutes, return the reaction to 4°C.


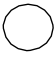




STEP	✓	Purify DNA	Notes
1		Add <b>1X</b> volume of AMPure PB beads to the End-Repair reaction.	
2		Mix the bead/DNA solution by tapping the tube.	
3		Allow the DNA to bind by letting it sit at room temperature for 10 minutes.	
4		Spin down the LoBind tube (for 1 second) to collect beads.	
5		Place the LoBind tube in a magnetic bead rack to collect the beads to the side of the tube.	
6		Slowly pipette off cleared supernatant and save (in another tube). Avoid disturbing the bead pellet.	
7		Wash beads with freshly prepared 70% ethanol.	
8		Repeat <a href="#">step 8</a> above.	
9		Remove residual 70% ethanol and dry the bead pellet. <ul style="list-style-type: none"> <li>– Remove the LoBind tube from the magnetic bead rack and spin to pellet beads. Both the beads and any residual 70% ethanol will be at the bottom of the tube.</li> <li>– Place the LoBind tube back on magnetic bead rack.</li> <li>– Pipette off any remaining 70% ethanol.</li> </ul>	
10		Check for any remaining droplets in the tube. If droplets are present, repeat <a href="#">step 10</a> .	
11		Remove the LoBind tube from the magnetic bead rack and allow beads to air-dry (with LoBind tube caps open) for 30 to 60 seconds.	
12		Elute the DNA off the beads in <b>30 µL</b> Elution Buffer. Mix by gently tapping the LoBind tube until homogenous, then let stand at room temperature for 2 minutes.	
13		Optional: Verify your DNA amount and concentration using a Nanodrop or Qubit quantitation platform, as appropriate.	
14		Optional: Perform qualitative and quantitative analysis using a Bioanalyzer system instrument with the DNA 12000 Kit. Note that typical yield at this point of the process (following End-Repair and one <b>1X</b> AMPure PB bead purification) is approximately <b>80-100%</b> of the total starting material.	
15		The End-Repaired DNA can be stored overnight at 4°C or at -20°C for longer duration.	

## Prepare Blunt-Ligation Reaction



Use the following table to prepare your blunt-ligation reaction:

1. In a LoBind microcentrifuge LoBind tube (on ice), add the following reagents in the order shown. Note that you can add water to achieve the desired DNA volume. If preparing a Master Mix, ensure that the adapter is NOT mixed with the ligase prior to introduction of the inserts. Add the adapter to the well with the DNA. All other components, including the ligase, should be added to the Master Mix.

Reagent	Tube Cap Color	Stock Conc.	Volume	Final Conc.	✓	Notes
DNA (End Repaired)			29.0 µL to 30.0 µL			
Blunt Adapter (20 µM)		20 µM	1.0 µL	0.5 µM		
Mix before proceeding						
Template Prep Buffer		10 X	4.0 µL	1X		
ATP low		1 mM	2.0 µL	0.05 mM		
Mix before proceeding						
Ligase		30 U/µL	1.0 µL	0.75 U/µL		
H <sub>2</sub> O			___ µL to adjust to 40.0 µL			
<b>Total Volume</b>			<b>40.0 µL</b>			

2. Mix the reaction well by gentle mixing.
3. Spin down contents of LoBind tube with a quick spin in a microfuge.
4. Incubate at 25°C for 15 minutes. At this point, the ligation can be extended up to 24 hours or cooled to 4°C (for storage of up to 24 hours).
5. Incubate at 65°C for 10 minutes to inactivate the ligase, then return the reaction to 4°C. You must proceed with adding exonucleases after this step.

## Exo III and Exo VII Treatment

Reagent	Tube Cap Color	Stock Conc.	✓	Volume
Ligated DNA				40 µL
Mix reaction well by pipetting				
ExoIII		100.0 U/µL		1.0 µL
ExoVII		10.0 U/µL		1.0 µL
<b>Total Volume</b>				<b>42 µL</b>

1. Spin down contents of LoBind tube with a quick spin in a microfuge.
2. Incubate at 37°C for 1 hour, then return the reaction to 4°C. You must proceed with purification after this step.

## Purify SMRTbell Templates

STEP	Purify SMRTbell Templates	Notes
1	Add <b>1X</b> volume of AMPure PB beads to the exonuclease-treated reaction.	
2	Mix the bead/DNA solution by tapping the tube.	
4	Allow the DNA to bind to beads by letting the sample sit at room temperature for 10 minutes.	
4	Spin down the LoBind tube (for 1 second) to collect beads.	
5	Place the LoBind tube in a magnetic bead rack to collect the beads to the side of the tube.	
6	Slowly pipette off cleared supernatant and save (in another tube). Avoid disturbing the bead pellet.	
7	Wash beads with freshly prepared 70% ethanol.	
8	Repeat <a href="#">step 8</a> above.	
9	Remove residual 70% ethanol and dry the bead pellet. <ul style="list-style-type: none"> <li>– Remove the LoBind tube from the magnetic bead rack and spin to pellet beads. Both the beads and any residual 70% ethanol will be at the bottom of the tube.</li> <li>– Place the LoBind tube back on the magnetic bead rack.</li> <li>– Pipette off any remaining 70% ethanol.</li> </ul>	
10	Check for any remaining droplets in the tube. If droplets are present, repeat <a href="#">step 10</a> .	
11	Remove the LoBind tube from the magnetic bead rack and allow beads to air-dry (with LoBind tube caps open) for 60 seconds.	
12	Elute the DNA off the beads in <b>50 µL</b> of Elution Buffer. Mix thoroughly by gently tapping the LoBind tube and let sit at room temperature for 2 minutes.	
13	The eluted DNA in <b>50 µL</b> Elution Buffer should be taken into the second <b>1X</b> AMPure PB bead purification step.	

STEP	Second Purification	Notes
1	Add <b>1X</b> volume of AMPure PB beads to the <b>50 µL</b> of eluted DNA.	
2	Mix the bead/DNA solution by tapping the tube.	
4	Allow the DNA to bind to beads by letting the sample sit at room temperature for 10 minutes.	
4	Spin down the LoBind tube (for 1 second) to collect beads.	
5	Place the LoBind tube in a magnetic bead rack to collect the beads to the side of the tube.	
6	Slowly pipette off cleared supernatant and save (in another tube). Avoid disturbing the bead pellet.	
7	Wash beads with freshly prepared 70% ethanol.	
8	Repeat <b>step 8</b> above.	
9	Remove residual 70% ethanol and dry the bead pellet. <ul style="list-style-type: none"> <li>– Remove the LoBind tube from the magnetic bead rack and spin to pellet beads. Both the beads and any residual 70% ethanol will be at the bottom of the tube.</li> <li>– Place the LoBind tube back on the magnetic bead rack.</li> <li>– Pipette off any remaining 70% ethanol.</li> </ul>	
10	Check for any remaining droplets in the tube. If droplets are present, repeat <b>step 10</b> .	
11	Remove the LoBind tube from the magnetic bead rack and allow beads to air-dry (with LoBind tube caps open) for 60 seconds.	
12	Elute the DNA off the beads in <b>50 µL</b> of Elution Buffer. Mix thoroughly by gently tapping the LoBind tube and let sit at room temperature for 2 minutes.	
13	Elute the off the beads in <b>10 µL</b> . Mix thoroughly.	
14	Determine concentration using a Qubit device or similar quantification assay.	
15	Run 1 µL of sample on Agilent DNA 12000 chip according to manufacturer's instructions. Perform qualitative analysis using a Bioanalyzer instrument with the DNA 12000 kit. Refer to Agilent Technologies' guides for specific information.	

## DNA Control Complex Dilution

You must have the PacBio Control Complex for this step. Dilute the Control Complex according to the volumes and instructions specified in Sample Setup.

## Anneal and Bind SMRTbell Templates

Follow the instructions in Sample Setup to anneal and bind your library.

Revision History (Description)	Version	Date
Initial release. Converted from "Unsupported Protocol" with updates to the cDNA preparation section.	01	June 2018

For Research Use Only. Not for use in diagnostic procedures. © Copyright 2013 - 2018, Pacific Biosciences of California, Inc. All rights reserved. Information in this document is subject to change without notice. Pacific Biosciences assumes no responsibility for any errors or omissions in this document. Certain notices, terms, conditions and/or use restrictions may pertain to your use of Pacific Biosciences products and/or third party products. Please refer to the applicable Pacific Biosciences Terms and Conditions of Sale and to the applicable license terms at <https://www.pacb.com/legal-and-trademarks/terms-and-conditions-of-sale/> Pacific Biosciences, the Pacific Biosciences logo, PacBio, SMRT, SMRTbell, Iso-Seq, and Sequel are trademarks of Pacific Biosciences in the United States and/or certain other countries. All other trademarks are the sole property of their respective owners.